



EZ-fill® Vial Platform



Clean. Sterile. Ready to Be Filled.

Stevanato Group is the largest manufacturer of sterilized vials. EZ-fill® ready-to-use vials are washed and sterilized, making them suitable for processing on both new and existing fill & finish lines, helping streamline pharmaceutical processes.




Features and Benefits

- **Standardize** - Market standard secondary packaging powered by SG across multiple / global vial suppliers
- **Increase flexibility** - EZ-fill® supports flexible fill & finish equipment for different packaging solutions and formats
- **Faster time to market** - Accelerate development from clinical to commercial with reliable ready-to-fill vial performance
- **Maintain** the highest quality - High quality manufacturing ensures consistent vials performances, allowing pharma companies de-risk operations and ease their manufacturing process
- **Reduce TCO** - Reduced total cost of ownership (TCO) driven by improved efficiency of the aseptic filling, reduced downtime and investments

Preferred Market Choice

- Vaccines
- Small Molecules
- Diluents
- Biologics
- ADCs
- Next Generation Lyo Drugs
- Animal Health

A PROVEN PROCESS

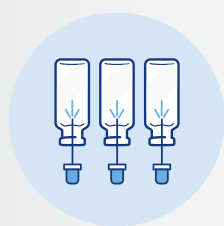

Washing



Depyrogenation



**Nesting & Tub
Insertion,
Tray Loading**



Sealing & Bagging



Sterilization

VIAL PORTFOLIO

	FORMAT		FILLING VOLUME	CRIMP FINISH	PCS
Nest & Tub* 	2R	2" TUB	3ml	13mm	120
	4R		4ml	13mm	120
	6R*		8ml	20mm	48
	8R		8ml	20mm	48
	10R*		10ml	20mm	48
	2R	3" TUB	3ml	13mm	100
	6R*		8ml	20mm	48
	10R*		10ml	20mm	48
	15R		15ml	20mm	24
	20R		20ml	20mm	24
	25R		25ml	20mm	24
	30R		30ml	20mm	24
	50R		50ml	20mm	16
Tray 	2R		3ml	13mm	228
	4R		4ml	13mm	228
	6R		8ml	20mm	96
	8R		8ml	20mm	96
	10R		10ml	20mm	96
	20R		20ml	20mm	60
	25R		25ml	20mm	60
	30R		30ml	20mm	60
	50R		50ml	20mm	28

ADDITIONAL FEATURES AVAILABLE

Baked-on siliconization: Prevents fogging and enables high recovery after reconstitution of lyophilized drugs
LDP: Low Delamination Propensity vials with optimized inner surface morphology for enhanced chemical durability
Amber glass: Shields light-sensitive drug products from degradation

VIAL TRAY CONFIGURATION

- Standard
- Upside-down


* The nest layouts differ between the two tubs

PACKAGING


Vials are delivered in nest and tub and tray configurations, in protective bag compatible with EtO sterilization




Protective Bag
Available in single- or double-bag configurations, depending on the selected de-bagging process

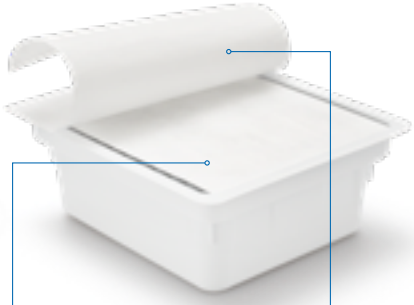

Liner
Provides protection for nested containers during lid opening

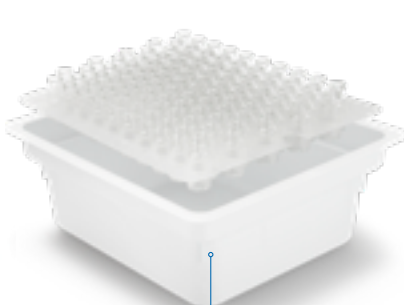

Sealing Lid
Sealed to the tub and designed for an easy and smooth peel-off

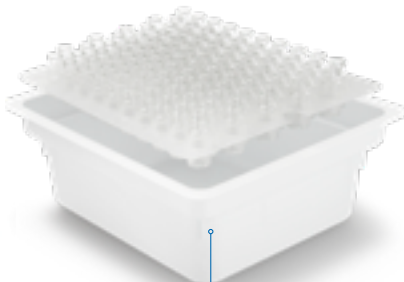

Nest & Tub
Designed to support different filling technologies

Headquarters

Via Molinella, 17
35017 Piombino Dese
Padova, Italy

**Global presence with operating units
and commercial offices in 10 countries.**

Visit our website to contact us.



stevanatogroup.com

SG, Stevanato, EZ-fill, EZ-fill Smart, Alba, Nexa, Fina, Alina, Vertiva are registered trademarks of Stevanato Group and/or its affiliates in one or more jurisdictions.

EZ-FILL_VIALS_202609